

Portable Meters

Denver Instrument Portable Meters go where you go!

Lightweight, ergonomically designed, waterproof meters with microprocessor accuracy

These models in Denver's Portable pH meter lineup share a rugged, ergonomic design, with a sealed waterproof case that's buoyant and fully immersible to 1 meter.

All models feature BNC and pin connectors, plus a waterproof Twist-Lock connector for special Denver Instrument electrodes.



Operation is exceptionally. The case fits comfortably in either hand, with soft rubber surfaces for a firm grip.

All measure pH with superb accuracy, with auto or manual temperature correction. And each one is ready to be toted most anywhere you go, providing up to 240 hours of use for Basic and about 40 hours for Advanced models on a set of four AA batteries. A selectable auto-off function saves batteries.

Three basic models offer accurate pH plus simple operation

- AP-5 for pH only
- AP-10 for pH, mV and relative mV readings

To operate, just press "Standardize" and the meter automatically recognizes any of 22 different buffers and checks electrode quality. You can calibrate with up to three buffers. Simple prompts and messages make accuracy easy.

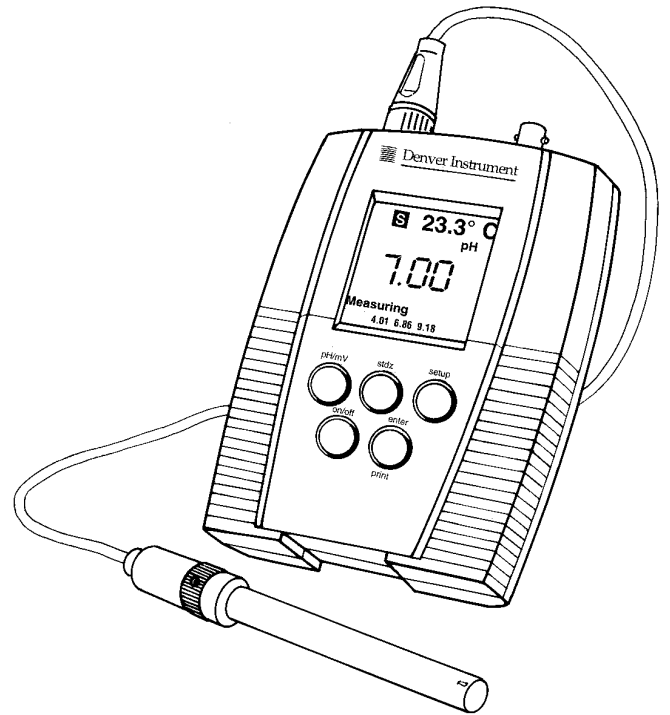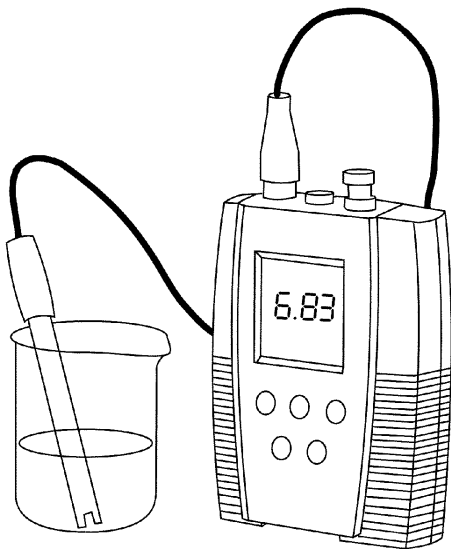

Specifications

pH	Range	-1.99 to +19.99
	Resolution	0.01
	Accuracy	±0.01
mV	Range	-999 to +999 (AP10 and AP15 only)
	Resolution	1
	Accuracy	±1mV
Temp	Range	-5 to +105.0°C
	Resolution	0.1°C
	Accuracy	±0.4°C
Buffer Recognition	Auto-recognition, 22 standard buffers	
Calibration	Up to 3 pH buffers	
Dimensions	4-7/8W x 6-5/8L x 1-3/4"H (12.5 x 17 x 4.5cm)	
Power Requirements	4 AA Batteries (included)	

Ordering Information

Meter with kit includes: Meter, plastic-body pH/ATC "3-in-1" electrode, softcase, 4 AA batteries and operation manual.
Meter only includes: Meter, 4 AA batteries, soft case and operation manual.

	AP5	AP10
Meter Kit	02-226-95	02-226-96
Meter Only	02-226-100	02-226-101

 **Denver Instrument Company**

Precision laboratory instruments since 1880
6542 Fig Street • Arvada, Colorado 80004
303-431-7255 • 1-800-321-1135 • Fax 303-423-4831
www.denverinstrument.com